## Instructions for Use

Bulgarian (BG)

Croatian (HR)

Czech (CS)

Danish (DA)

English-UK (GB)

Estonian (ET)

French (FR)

German (DE)

Greek (EL)

Hungarian (HU)

Italian (IT)

Lithuanian (LT)

Macedonian (MK)

Norwegian (NO)

Polish (PL)

Portuguese (PT)

Romanian (RO)

Russian (RU)

Slovak (SK)

Spanish (ES)

Swedish (SV)

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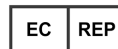
EDX Ehrlichia Run Control Summary	
Catalogue number	ERLRC
Analyte	<i>Ehrlichia chaffeensis</i>
Matrix	Synthetic Matrix
Preservative	ProClin® 300
Storage	-20 °C or below
Fill volume	1.0 mL
Number of vials per kit	5 vials
Number of uses	Multiple use
Shelf life	18 months from date of manufacture
Precautions	Biological risks
Regulatory status	For in vitro diagnostic use

**INTENDED USE:**

EDX Ehrlichia Run Control is an external quality control intended to be used with *Ehrlichia chaffeensis* assays to monitor the performance of that given assay. The control is used to monitor the presence of *Ehrlichia chaffeensis* DNA. Routine use of EDX Ehrlichia Run Control allows laboratories to evaluate day-to-day and lot-to-lot variation of their molecular assay and test for operator proficiency.

**PRODUCT DESCRIPTION:**

EDX Ehrlichia Run Control contains whole, intact organisms. The control is intended to test the entire process of a molecular assay including extraction, detection, and amplification.

EDX Ehrlichia Run Control is internally value assigned using Bio-Rad Digital Droplet PCR (ddPCR). EDX Ehrlichia Run Control is formulated in Synthetic Matrix that contains preservatives.

**PROCEDURE:**

EDX Ehrlichia Run Control requires an extraction step. The control must be treated in a similar manner to other tested samples. This allows the operator to assess their extraction technique.

EDX Ehrlichia Run Control is an unassayed control and does not have assigned values. The control is not a substitute for the manufacturer's kit controls provided with the assay. EDX Ehrlichia Run Control is to be tested according to the assay manufacturer's or testing laboratory's instructions and recommendations. Expected results using EDX Ehrlichia Run Control must be established by the end user for their specific assay.

**STORAGE AND HANDLING:**

EDX Ehrlichia Run Control should be stored at -20 °C or below.

Thaw EDX Ehrlichia Run Control at room temperature and vortex briefly prior to use.

Once thawed, EDX Ehrlichia Run Control is stable for 24 hours when stored at 2-8 °C. Product can be frozen and reused a maximum of four (4) times. After the fourth freeze-thaw, any remaining materials must be disposed of. Do not dilute.

Do not use EDX Ehrlichia Run Control beyond the expiration date.

**LIMITATIONS:**

For in vitro diagnostic use.

**WARNINGS AND PRECAUTIONS:**

EDX Ehrlichia Run Control should be considered biohazardous. Universal precautions and proper disposal should be practised<sup>1</sup>. Do not pipette by mouth. Do not smoke, eat, or drink in areas where specimens are handled.

Discard product if packaging is damaged or leaking. Disinfect liquids, materials, or spills with a 0.5% sodium hypochlorite solution. Dispose of all materials and liquids used in the procedure as if they contain pathogenic agents.

**REFERENCES:**

<sup>1</sup> Siegel JD, Rhinehart E, Jackson M, Chiarello L, and the Healthcare Infection Control Practices Advisory Committee, 2007  
Guideline for Isolation Precautions: Preventing transmission

BIO-RAD and DROPLET DIGITAL PCR are trademarks of Bio-Rad Laboratories, Inc. in certain jurisdictions.

**SYMBOL REFERENCES:**

Catalogue number

Lot number

For in vitro diagnostic use

Biological risks

Expiration date

Upper limit of temperature

Manufacturer

Caution

Consult instructions for use


Positive Control

European Mark of Conformity

Authorised Representative

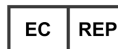Emergo Europe  
Prinsessegracht 20  
2514 AP The Hague  
The Netherlands

Exact Diagnostics LLC  
3400 Camp Bowie Blvd CBH-214  
Fort Worth, Texas 76107  
United States

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BIO-RAD och DROPLET DIGITAL PCR är varumärken som tillhör Bio-Rad Laboratories, Inc. i vissa länder.